- WATER FLOSSER • MUNDDUSCHE • ИРРИГАТОР ПОЛОСТИ РТА •
- IRRIGATEUR BUCCAL •

MODEL/MODELL
МОДЕЛЬ/ MODÈLE: VIP-008

VES electric

**EN INSTRUCTION MANUAL
DE BEDIENUNGSANLEITUNG
RU РУКОВОДСТВО ПО ЭКСПЛУАТАЦИИ
FR MODE D'EMPLOI**

WATER FLOSSER

MODEL: VIP-008

Dear customer!

Before using the appliance read the instruction manual carefully. Keep the instruction manual safe and pass it on if you pass the appliance on.

Technical characteristics:

Max. water pressure: 900 kPa / Power supply: Li-ion accumulator 3.7 V / Charger rating: 100-240 V 50/60 Hz 1.4 W / Frequency: 1700 pulses / min max. / 3 operation modes / Water tank capacity: up to 170 ml

Recommendations for the buyer

When buying the appliance, check it for the absence of mechanical damage of the construction, check the external design and packing, the functioning, the presence of all accessories in set. Note: operation checking - is the duty of the seller. Please check the warranty card and the accuracy of its filling (presence of the seller's seal). You must save the bill, instruction manual, warranty card until the end of the warranty period. In case of bill and warranty card loss the warranty time is counted from the production date.

PRECAUTIONS

- Use only the filtered water! Tap water contains salt, sand, lime, etc., which will clog the appliance. The clogged / dirty appliance or the appliance with poor cleaning voids the warrantee.

- Water temperature must not exceed 40°C!

• Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision. • the appliance is for home use only. Use the appliance for the direct purpose. • use the appliance strictly in accordance with your dentist recommendations. • use the appliance only with the supplied adapter and tips. • before the charging make sure that the voltage in your wall outlet is the same as on the charging adapter's rating label. • if you dropped the appliance / adapter and it has some visible defects, discontinue its use, bring to an authorized service center for checking. • do not drop or kick the appliance. • do not immerse the appliance / adapter or adapter's cord in water. • people that can not operate the appliance themselves, or people with low oral sensitivity should not use this appliance. • keep the tips far from children. The tips are small enough to be

swallowed and can become lodged in the throat. • instruct and look after the children or people with physical incapacibilities to use the appliance correctly. • do not let working appliance be unattended. • do not use any abrasive, chemical cleaners, polish for cleaning the appliance. • do not use the damaged appliance, damaged adapter or its cord. • do not let the adapter's cord hang over the edge of table or touch any hot surfaces. • repair the appliance only in the authorized service center. Do not repair the appliance by yourself! If the opening of the appliance is detected, the appliance's guarantee will be canceled. • while using the appliance, please, observe the safety precautions. • incorrect operation and improper handling can lead to malfunction of the appliance or injuries to the user.

• do not leave the appliance where it can drop in water. If the appliance fell in the water, discontinue the use. Bring it to an authorized service center for checking. • do not direct the water jet under your tongue, into the ear, nose or at any body part. Water pressure can cause serious harm to body parts. • fill the water tank with water or fluids (rinsing or medicinal fluids), recommended by your dentist and nothing else. The liquid should not leave sediment or be viscous, this may soil the appliance. • in case you are having oral treatment, or suffer with severe paradontosis consult the dentist before using the appliance. • do not insert foreign objects into the appliance.

• do not use bleaching or iodine-containing solvents. • it is necessary to pour water remnants out of water tank, then leave the appliance with the opened water tank cap until the tank is fully dry otherwise mold and bacteria can appear in the appliance and can affect your health, lead to unpleasant smell of the appliance and worsen appliance's work. • do not turn on the appliance with empty water tank. • don't use the appliance in the vicinity of hot surfaces, places where explosive gas, paint thinners, gasoline, paint or other flammable liquids are used or stored. • do not use this appliance if you have jewelries in the mouth; take them away before using the appliance. • do not use the appliance if you have open sore on tongue / in the mouth.


• give a high attention when the appliance is being used near the children. Do not let children use or play with the appliance.

• do not use the appliance outdoors.

DESCRIPTION

1. tip / 2. body / 3. on/off button / 4. mode button
5. mode indicators / 6. charging indicator / 7. charging

base / 8. tip release button / 9. suction tube (don't detach the tube!)
 10. water inlet / 11. water tank / 12. water tank cap
 13. filter (don't detach the filter!)

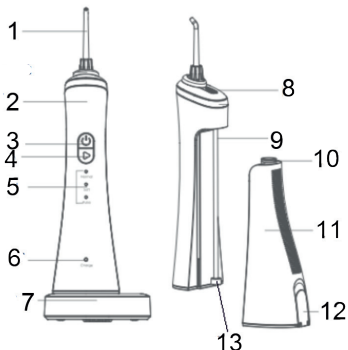Note:
 There may be some liquid remaining inside the appliance before use (tank, suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.


BEFORE FIRST USE

Remove the appliance from the package. Do not leave the packing (plastic bags, polystyrene etc.) within the reach of children since they could be dangerous. Keep the package for appliance storage.

After transportation in low temperatures leave the appliance in the room for 2 hours before use.

CHARGING

The manufacturer does not guarantee full charge of accumulator of the new appliance, that is why charge the appliance for about 24 hours before the use. The next charges will last 15-16 hours.

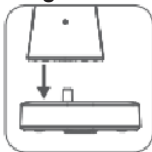Insert the charging base securely into the power outlet. Place the appliance on the charging base. The red light of Charge indicator will start blinking. When the Charge light goes blue light means the battery is fully charged.

Note: If the Charge indicator is red blinking during or after an irrigating, that means the

battery needs to be recharged.

MAINTENANCE

Note: If your physician recommended you to take antibiotics before having oral hygiene procedures, consult your dentist before using given appliance or any other appliance for oral

hygiene.

Water tank filling

Water temperature must not exceed 40°C! Use only the filtered water! Tap water contains salt, sand, lime, etc., which will clog the appliance. The clogged / dirty appliance or the appliance with poor cleaning voids the warrantee.

Open the water tank cap and fill water tank with warm water, close the lid. Water tank capacity is 170 ml. Do not overfill the tank!

Attention! Do not use bleaching or iodine-containing fluids.

Fill the water tank with water or fluids (rinsing or medicinal fluids), recommended by your

dentist and nothing else. The liquid should not leave sediment or be viscous, this may soil the appliance. If you used fluid, you should pour out its remnants

and clean the appliance. For this fill the water tank with filtered water and drain water through used tip, turning the appliance on. Make cleaning over the sink.

Tip installation and detachment

Install and detach tip with switched off appliance only. **To detach** the tip press the tip release button. There will be the click sound. Pull the nozzle upwards. **Insert** the tip into the tip inlet, press down until the click sound. Do not refit the tips while the appliance is working. Tips are 360° rotatable. Tips are the items of personal hygiene

and do not have warranty.


Functional buttons Use on/off button to switch the appliance on and off. Press mode button to select the mode you desire. During irrigating, you can change the working mode.

Normal: for everyday use. / Soft: for gentle cleaning of sensitive teeth and gums. / Pulse: for massaging gums to promote blood circulation and strengthen the teeth.

Location in the mouth It is not allowed to use the appliance at the angle of 180° or turning it upside down. Use the appliance at the angle of 45°- 90°.

Locate the tip in the mouth. Direct the tip end at the teeth. Switch the appliance on. Set the required water level. Direct water jet, turning the appliance or tip. Close your lips to avoid water splashing, but do not close your mouth allowing

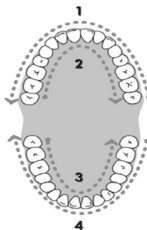water freely flow away. For better result please start with side teeth and molars, moving it to fore teeth. Move tip from tooth to tooth making pauses between moves. Go on cleaning until teeth surrounding zones and zones between teeth are clean. **The appliance will shut off automatically after operating for 2 minutes. Please press the on/off button again if need to continue cleaning.**

TIPS There are 3 standard tips in set to be used for teeth cleaning. The tips can be stored on the charging base.

CLEANING 1. Switch the appliance off. 2. Open the water tank cap. Pour water remnants out of water tank. Rinse water tank and let dry with the cap open. You can use 40% medicinal spirit solution to disinfect the appliance. Pour the solution to the tank and drain through the used tip, turning the appliance on. After disinfection, drain the clean the water. Dry the tank after each use. 3. Wipe the body with wet cloth, wipe dry. Do not use chemical, abrasive substances or polish to clean the appliance, do not use wire brushes. Don't wash the charging base!

STORAGE Store the appliance empty in clean place, far from children. While transportation of appliance it should be dry and empty.

ATTENTION! *Contact your dentist, if: you doubt the safety of appliance use; you have intense gum bleeding after the appliance use, or gum bleeding still exists 2 weeks since the beginning of the appliance use.*

TROUBLE SHOOTING

Trouble	Purpose	Trouble shooting
The appliance does not work.	You've just purchased the appliance or it has been left unused for more than 3 months.	Charge the appliance for 24 hours.
The appliance works only for a few minutes even after charging.	Charging time was insufficient.	Charge the appliance for 24 hours.

Low water pressure	The dirt stuck in the tip. The tip is deformed. The filter is blocked.	Clean the tip. Replace the tip. Clean the filter.
Water does not come out.	The tank is empty. The appliance is tilted too much while in use.	Fill up the tank. Use the appliance in a vertical state.

NOTE! Manufacturer reserves the right to change technical characteristics of the device, with the aim of improving its quality, also change the parts and color of the appliance.

Disposal

Dispose the appliance in an environmentally friendly manner. It should not be disposed with the normal household waste. Please dispose it at a recycling centre for used electrical and electronic devices. You can obtain more information from your local authorities.

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